

# Phonak Audeo™ Fit

## Product information

Our latest innovation in the Lumity line, Audeo Fit has a heart rate sensor that can track the body's health data which can be viewed with the myPhonak app. By tracking and synchronizing health data such as calories burned, heart rate and activity levels, these hearing aids can promote healthy habits and create a greater sense of well-being. Suitable for mild to severe hearing losses and powered by the innovative technology of the Phonak Lumity platform, Phonak Audeo Fit hearing aids are the solution for those clients that want to engage with their health data.





# Product description



Audéo L-R Fit  
Audéo L-R Fit Trial  
L90/L70

Multi-function button		•
Bluetooth®		•
Hydrophobic coating		•
IP rating		IP68 <sup>1</sup>
Dimensions	L x W x D	30.4 x 12.6 x 8.1mm 1.2 x 0.5 x 0.32
	Weight	0.25g (0.079 oz)
Charging times <sup>2</sup>	0% to 100%	3 h
	0% to 80%	1.5 h
	30 min. fitting	10 min.
<b>M Receiver</b>		
Maximum Power Output (dB SPL)		114
Max. gain (dB)	2 cc coupler	51
Frequency range (Hz)	2 cc coupler	>100 - <8000
<b>P Receiver</b>		
Maximum Power Output (dB SPL)		122
Max. gain (dB)	2 cc coupler	59
Frequency range (Hz)	2 cc coupler	>100 - <6300

## Housing colors



**P1**  
Sand Beige



**P4**  
Chestnut



**P6**  
Silver Gray

<sup>1</sup> IP68 indicates that the hearing aid is water resistant and dust tight. It survived continuous immersion in 1 meter of water for 60 minutes and 8 hours in a dust chamber as per the IEC60529 standard.  
<sup>2</sup> Actual charging time depends on the remaining battery life, but it will not exceed 3 hours.



# Performance levels

AutoSense OS™ 5.0	Premium	Advanced
Comfort in Echo	•	
Speech in Car	•	
Speech in Loud Noise	•	
Music	•	•
Comfort in Noise	•	•
Speech in Noise	•	•
Calm Situation	•	•
Media speech + mic	•	•
Media music + mic	•	•
<b>Additional programs</b>		
Max. additional programs	4	4
Comfort in Echo	•	
Speech in Loud Noise	•	•
Speech in 360°	•	•
Speech in Noise	•	•
Calm Situation	•	•
Comfort in Noise	•	•
Music	•	•
Acoustic phone	•	•
Custom program	•	•
<b>Streaming programs</b>		
RogerDirect + mic	•	•
PartnerMic + mic	•	•
Phone call + mic	•	•

Features	Premium	Advanced
SpeechSensor	•	
Speech Enhancer	•	
StereoZoom 2.0	•	•
Dynamic Noise Cancellation	•	•
Motion Sensor Hearing	•	•
RogerDirect™	•	•
Roger and Directional streaming	•	•
EchoBlock	•	
Tap control		
- Siri/Google Assistant	•	•
- Pause/resume streaming	•	•
- Accept/end phone call	•	•
SoundRelax	•	•
DuoPhone	•	•
Real Ear Sound	•	•
UltraZoom	•	•
Health functionalities	•	•
WindBlock	•	•
User Preference Tuning	•	•
Tinnitus Balance	•	•
Auto Acclimatization	•	•
Environmental balance	•	•
SoundRecover2	•	•
WhistleBlock	•	•
NoiseBlock	•	•
QuickSync	•	•
Fine tuning channel	20	20

SmartSpeech™ Technology

\*availability country dependent

# Health features

Health feature progression	Audéo L-R	Audéo L-R Fit
Wearing time	•	•
Sound environments	•	•
Step counter	•	•
Distance walked or ran	•	•
Daily Goals	•	•
Sync. with Apple Health (new)	•	•
Live heart rate		•
Daily heart rate graph		•
Resting heart rate		•
Calories burned (new)		•

The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc.  
 IOS is a trademark of Cisco, Inc.  
 Android is a trademark of Google, Inc.

# Fitting

## Audéo L-R Fit Audéo L-R Fit Trial

Phonak Target	9.0 or higher
Noahlink Wireless	•

## Acoustic coupling\*

### Length

Available in four different lengths (0, 1, 2, 3), left and right

### Dome

Cap dome

Open dome Available in three sizes (S, M, L)

Vented dome clear Available in three sizes (S, M, L)

### cShell

Acryl


### Fitting ranges

### M Sensor Receiver


### P Sensor Receiver

Available in four different lengths (0, 1, 2, 3), left and right	•	•
Cap dome	•	
Open dome Available in three sizes (S, M, L)	•	
Vented dome clear Available in three sizes (S, M, L)	•	
Acryl	•	•

\* Only the 8 pin Sensor Receiver and Sensor cShell works with Audeo Fit.

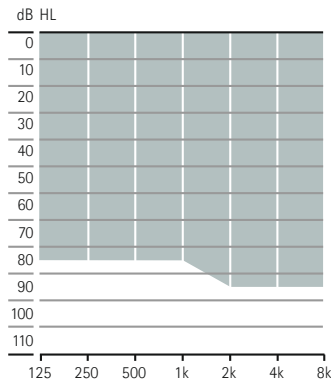

Mild to severe hearing loss.

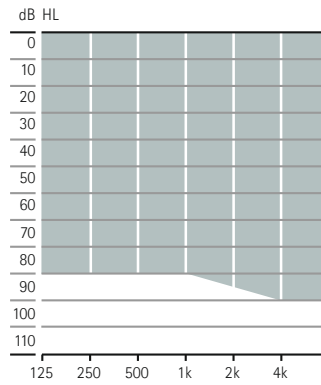

Mild to severe hearing loss.

# Phonak Audéo Fit charging option



## Phonak Charger Ease

The Phonak Charger Ease is a compact charger with magnets in the charging slots that make it easy to insert and remove Audéo L-R Fit. It uses a USB-C charging port that lights up when power is applied.

## Phonak Audéo Fit Bluetooth® versions

Direct calls/audio streaming

Bluetooth® 4.2 wireless technology and most older Bluetooth phones

The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc.

## Wireless communication portfolio

	Audéo L-R Fit Audéo L-R Fit Trial L90/L70
<b>Wireless accessory</b>	
TV Connector	•
Phonak PartnerMic	•
Phonak RemoteControl	•
<b>Roger</b>	
RogerDirect™	•
<b>Apps</b>	
myPhonak	•

For more information please visit [www.phonakpro.com/audeolumity](http://www.phonakpro.com/audeolumity)  
 For information on Roger microphones please visit [www.phonakpro.com/roger](http://www.phonakpro.com/roger)  
 For more information please visit [www.phonakpro.com/smartapps](http://www.phonakpro.com/smartapps)



# life is on


Since 1947, Phonak is dedicated to preserving life quality by opening new acoustic worlds. Already back then, in the days of our foundation, our company was driven by a formative conviction: We believe that well-hearing equates to well-being and thus is essential for living life to the fullest. In fact, the sense of hearing is directly linked to social, emotional, cognitive and physical well-being. Today as in future, we strive to offer the broadest portfolio of innovative hearing solutions. And, together with our hearing care professionals, we keep on focusing on what matters most: improving speech understanding, changing people's lives and having a positive effect on society as a whole.

[www.phonak.com](http://www.phonak.com)



027-0694-02/V1.00/2023-06/na © 2023 Sonova AG All rights reserved



 Sonova AG · Laubisrütistrasse 28  
CH-8712 Stäfa · Switzerland

